

Rehabilitation & Wellness
MINATO
MINATO MEDICAL SCIENCE CO., LTD.

Massage Robot

AQUATIZER

Dry Hydro Jet Massage Bed **QZ-240**

Japan
Good Design
Award




**GOOD
DESIGN**



※ illustrations are all image diagram.


The “four new nozzles” will provide a wide variety of treatment options.

Four nozzles

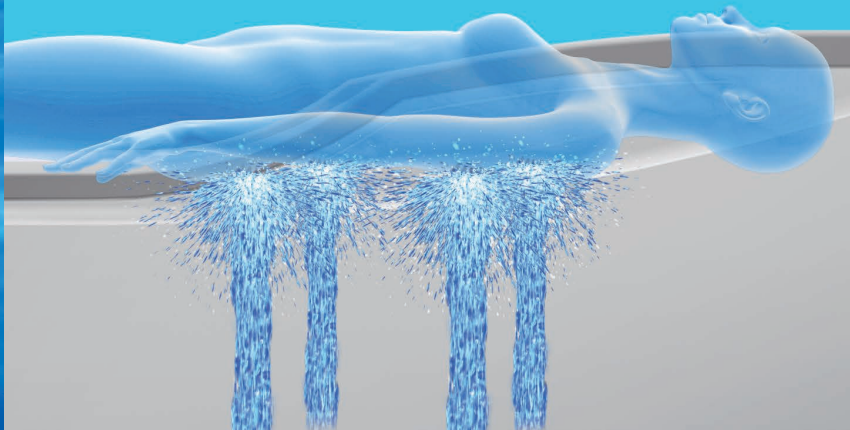




By alternating between four nozzles and two nozzles according to treatment areas, treatment patterns increase, such as "wide-ranging of treatment" and "treatment concentrated to a necessary area."

Two nozzles

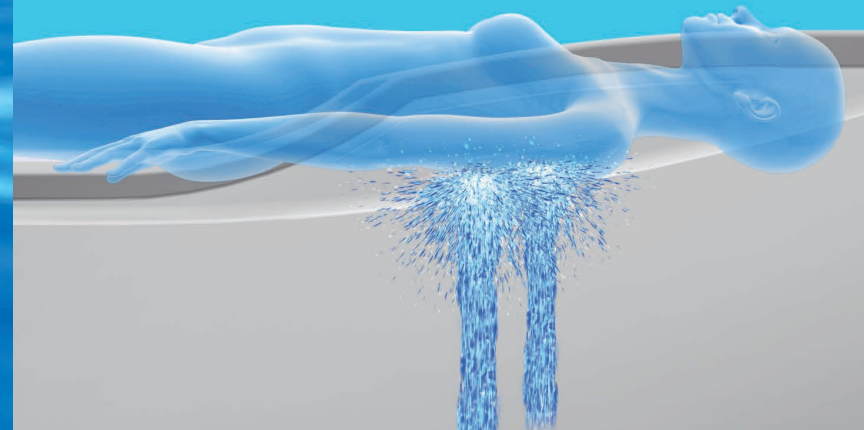

As is conventionally done, treatment can be performed using two nozzles only.

Various Massage Patterns

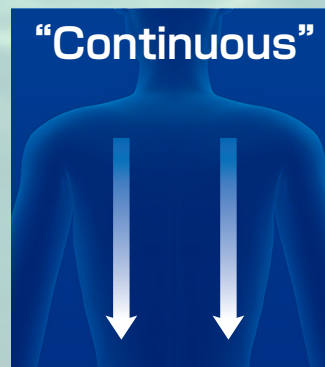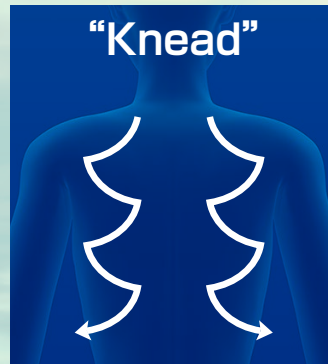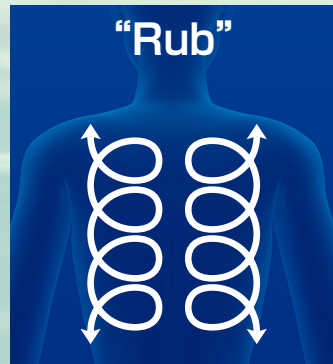

Reproducing
7 Types of
Skillful Manual
Massage.

Switching between the "Standard" and "Slow" modes in addition to the 7 massage patterns can provide more customized massage modes.

Foot Airbags, gently holding the legs



Unless the feet are held, the legs are open, and the massage area is deviated from the center of jet flow.



Because the legs are prevented from opening by softly holding the feet, the jet flow adequately catches the area from the calf to the thigh.

Foot Airbags, gently holding the legs



Thus far, as the feet were not held, the head made a swinging motion.



Holding the feet softly can stabilize the body and reduce the swinging motion of the head.




Newly Mounted Sensor for Automatic Head Position Detection

The newly mounted "height detection sensor" easily identifies the areas to be massaged.




Newly Mounted Sensor for Automatic Head Position Detection

When a patient lies on the bed, the sensor detects the head position, automatically measures the patient's height, and sets the area to be massaged according to the specified treatment site.

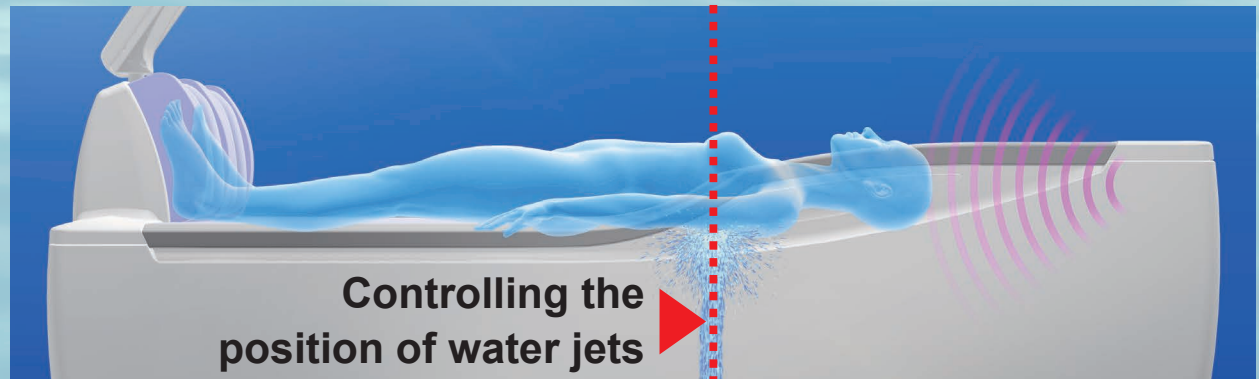


Excellent mat tension that directly carries the jet flow.


Softly stretched mat can carry the natural jet flow to the affected area as it is.

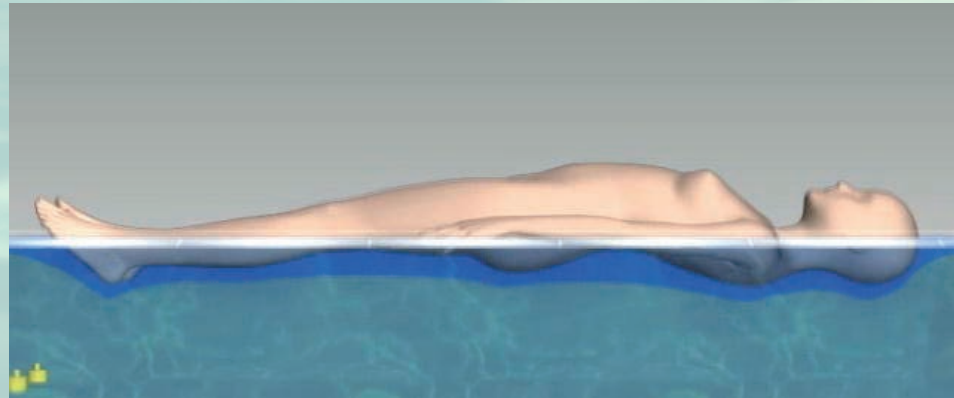



Excessively strong mat tension may produce an area to which the jet flow is not carried.

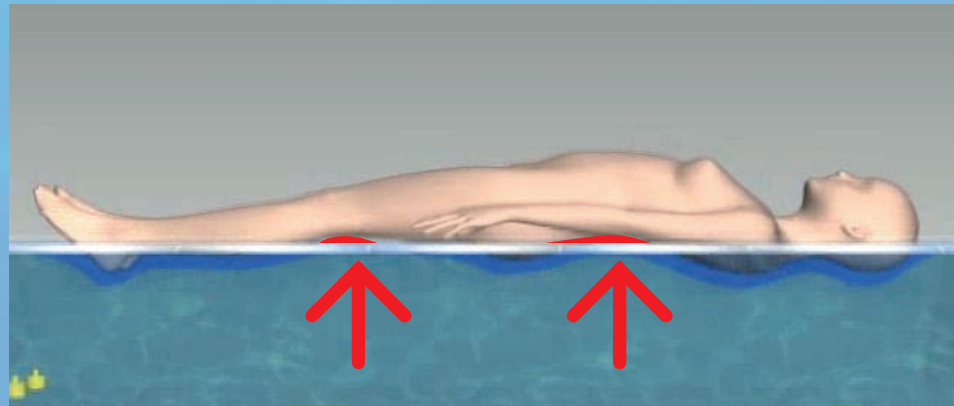


The "TURBO" mode provides a more powerful massage.

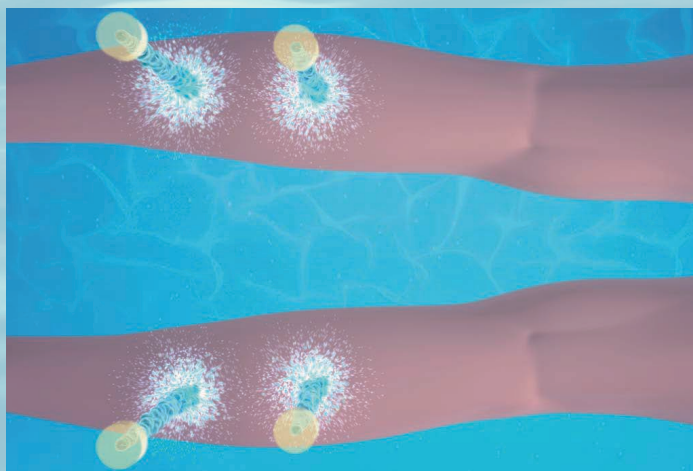


"NORMAL"

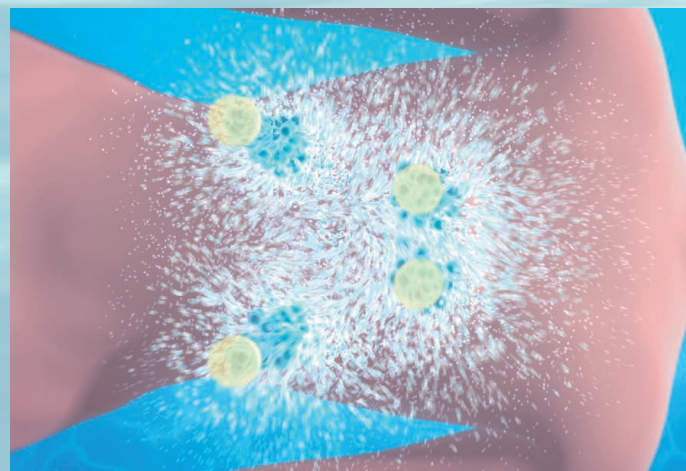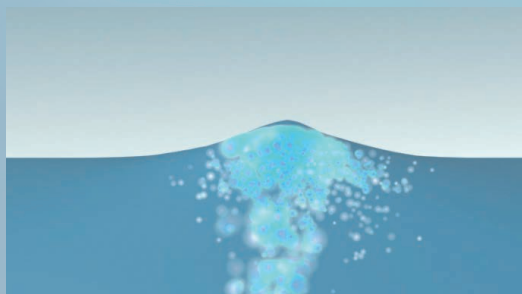

"TURBO"

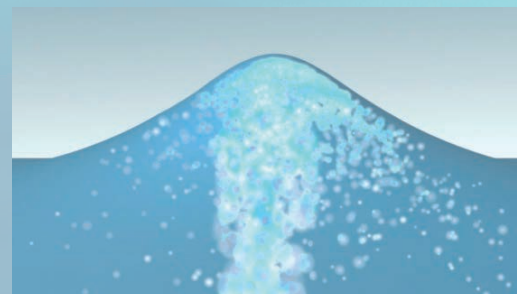


The lineup of the hooded type

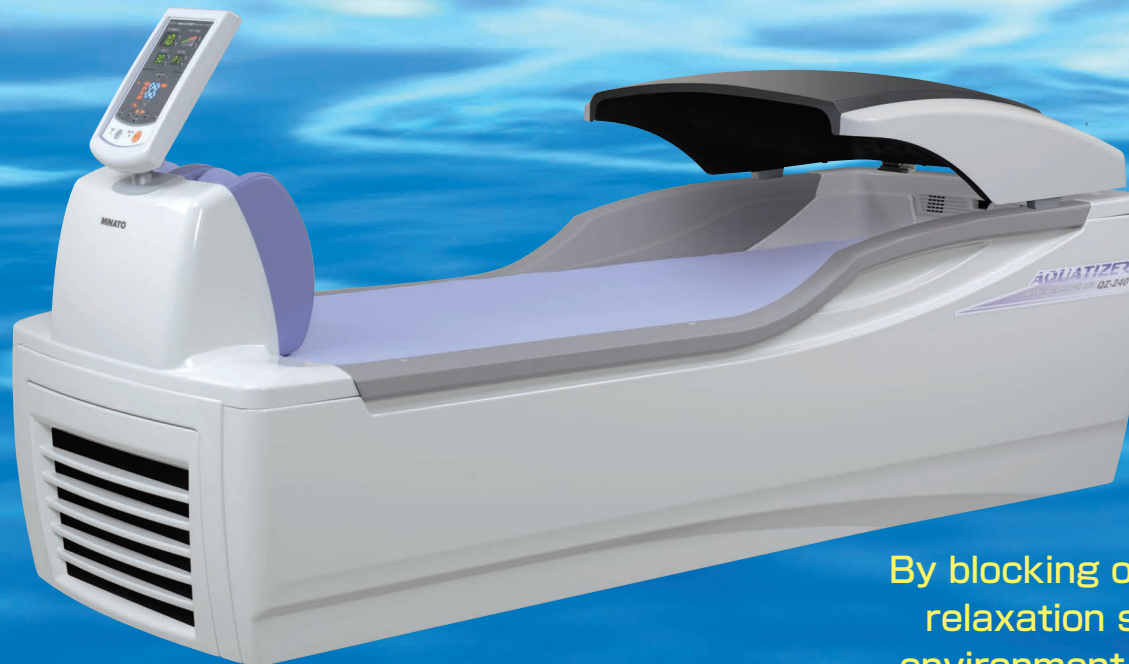

By blocking out the surroundings, a relaxation space is created in the environment and more comfortable treatment can be provided.

The lineup of color variation



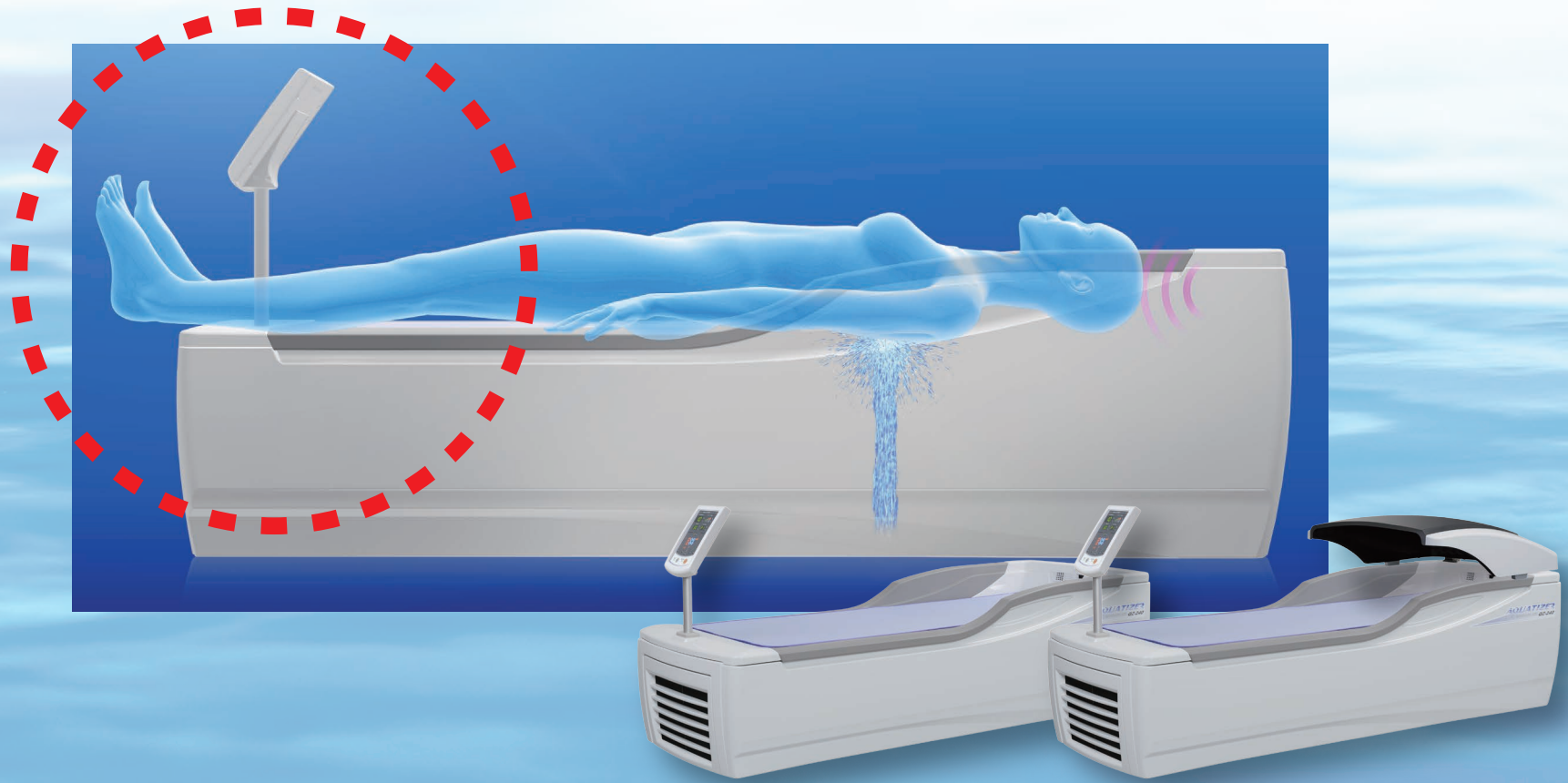


The type of Dry Hydro jet Massage Bed QZ-240 without the foot airbags is optionally available for the facilities where there are many persons of which height is over 190cm.




Next-generation Dry Hydro Jet Massage Bed

The **No. 1** market share in Japan

Changes in the dry hydro jet massage bed market in Japan

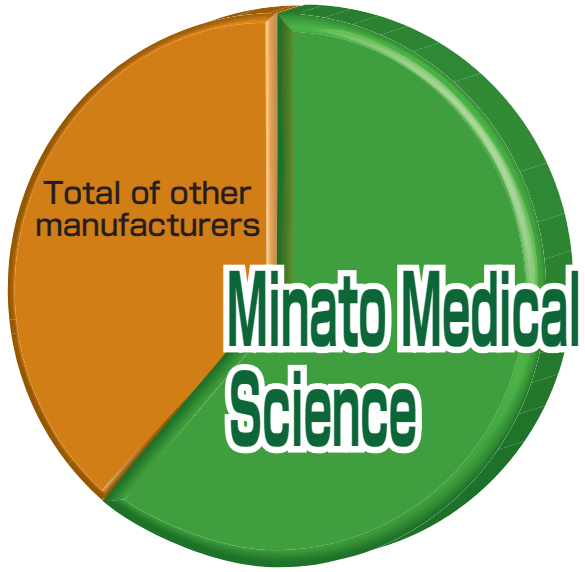

Market shares for
massage beds in Japan

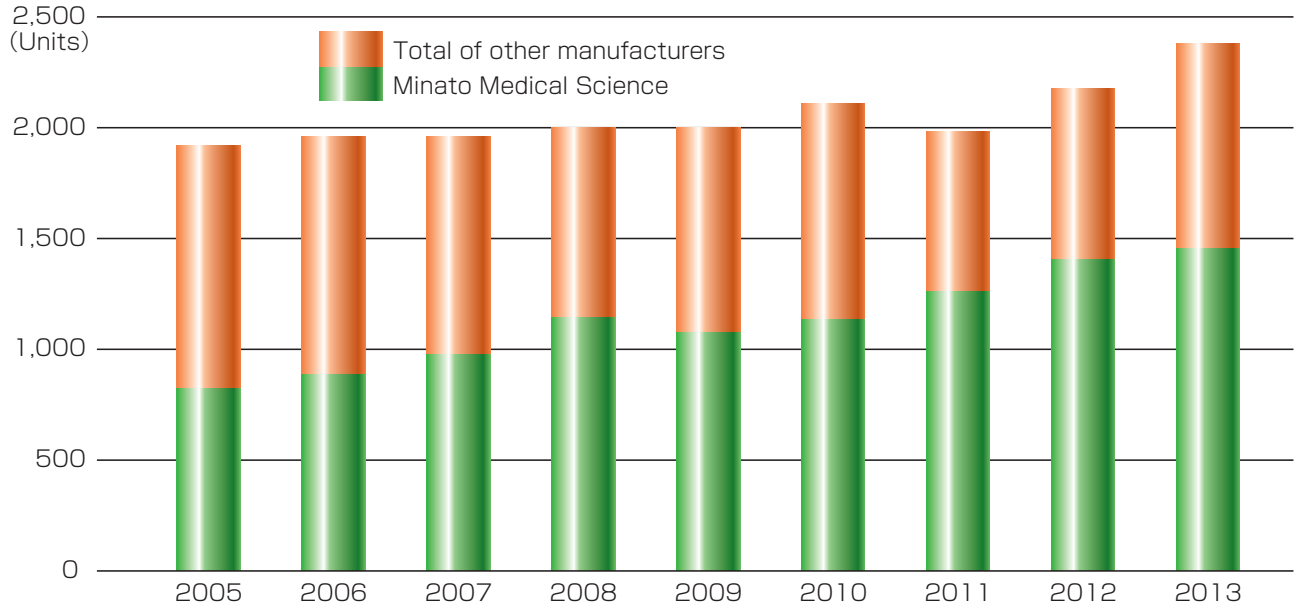

※Study by Yano Research Institute Ltd.